



ANALYTICAL BALANCES


PROFESSIONAL MEASURING





KERN Pictograms


 **Internal adjusting**
Quick setting up of the balance's accuracy with internal adjusting weight (motordriven)


 **Adjusting program CAL**
For quick setting up of the balance's accuracy. External adjusting weight required


 **EasyTouch**
Suitable for the connection, data transmission and control through PC or tablet


 **Memory**
Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.


 **Alibi memory**
Secure, electronic archiving of weighing results, complying with the 2014/31/EU standard.


 **KERN Universal Port (KUP)**
allows the connection of external KUP interface adapters, e.g. RS-232, RS-485, SB, Bluetooth, WIFI, Analogue, Ethernet etc. for the exchange of data and control commands, without installation effort


 **RS-232 Data interface**
To connect the balance to a printer, PC or network


 **RS-485 Data interface**
To connect the balance to a printer, PC or other peripherals. Suitable for data transfer over large distances. Network in bus topology is possible


 **USB Data interface**
To connect the balance to a printer, PC or other peripherals

 **Bluetooth® Data interface**
To transfer data from the balance to a printer, PC or other peripherals

 **WIFI Data interface**
To transfer data from the balance to a printer, PC or other peripherals

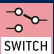
 **Control outputs**
(optocoupler, digital I/O)
To connect relays, signal lamps, valves, etc.


 **Analogue interface**
to connect a suitable peripheral device for analogue processing of the measurements

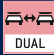
 **Interface for second balance**
For direct connection of a second balance


 **Network interface**
For connecting the scale to an Ethernet network


 **KERN Communication Protocol (KCP)**
It is a standardized interface command set for KERN balances and other instruments, which allows retrieving and controlling all relevant parameters and functions of the device. KERN devices featuring KCP are thus easily integrated with computers, industrial controllers and other digital systems


 **GLP/ISO log intern**
The balance displays weight, date and time, independent of a printer connection


 **GLP/ISO log Printer**
With weight, date and time. Only with KERN printers.


 **Piece counting**
Reference quantities selectable. Display can be switched from piece to weight


 **Recipe level A**
The weights of the recipe ingredients can be added together and the total weight of the recipe can be printed out


 **Recipe level B**
Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display


 **Totalising level A**
The weights of similar items can be added together and the total can be printed out


 **Percentage determination**
Determining the deviation in % from the target value (100 %)


 **Weighing units**
Can be switched to e.g. nonmetric units. See balance model. Please refer to KERN's website for more details

 **Weighing with tolerance range (Checkweighing)**
Upper and lower limiting can be programmed individually, e.g. for sorting and dosing. The process is supported by an audible or visual signal, see the relevant model


 **Hold function**
(Animal weighing program)
When the weighing conditions are unstable, a stable weight is calculated as an average value


 **Protection against dust and water splashes IPxx**
The type of protection is shown in the pictogram


 **Suspended weighing**
Load support with hook on the underside of the balance


 **Battery operation**
Ready for battery operation. The battery type is specified for each device


 **Rechargeable battery pack**
Rechargeable set


 **Universal plug-in power supply**
with universal input and optional input socket adapters for
A) EU, CH, GB
B) EU, CH, GB, US
C) EU, CH, GB, US, AUS


 **Plug-in power supply**
230V/50Hz in standard version for EU, CH. On request GB, USA or AUS version available


 **Integrated power supply unit**
Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, USA or AUS on request


 **Weighing principle Strain gauges**
Electrical resistor on an elastic deforming body


 **Weighing principle Tuning fork**
A resonating body is electromagnetically excited, causing it to oscillate


 **Weighing principle Electromagnetic force compensation**
Coil inside a permanent magnet. For the most accurate weighings


 **Weighing principle Single cell technology**
Advanced version of the force compensation principle with the highest level of precision

 **Conformity Assessment**
The time required for conformity assessment is specified in the pictogram

 **DAkkS calibration possible (DKD)**
The time required for DAkkS calibration is shown in days in the pictogram

 **Factory calibration (ISO)**
The time required for Factory calibration is shown in days in the pictogram

 **Package shipment**
The time required for internal shipping preparations is shown in days in the pictogram

 **Pallet shipment**
The time required for internal shipping preparations is shown in days in the pictogram

KERN Models A-Z

440	30
572	34
A	
ABP	50-51
ABP-A	52
ABT-NM	48
ACS/ACJ	46-47
ADB/ADJ	43
ALS-A/ALJ-A	44-45
B	
BFB	124-125
BFN	126
BIC	121
BID	122-123
C	
CB Q1 · CB Q2 · CB P1	160
CCA	94-95
CCS	96-97
CDS	93
CE Hx	138
CFS	90
CH	169
CIB	87
CJ P · CJ X	163
CKE	91
CM	10
CO Y1 · CO Y2 · CO Y5	162
CP P4 · CP Y4 · CP P1 · CP Y1 · CP P3	157
CP P2 · CP P8 · CP P7 · CP P9	158
CR Q1 · CR P1 · CR Y1	159
CT Q1 · CT P1 · CT P2	161
CPB	89
CXB/CXB-NM	88
D	
DAB	54
DBS	55
DE	106-107
DLB	56
DS	115
E	
ECB-N/ECE-N	59
EFS	13
EG-N/EW-N	37
EHA	17
EMB	14
EMB-V	15
EMS	16
EOB	104
EOC	108-109
EOE	103
EOS	105
EWJ	31

F	
FCB	61
FCE-N	60
FCF	69
FES/FEJ	41
FFN	74
FGE	68
FKB	62-63
FOB-LM	70
FOB-S · FOB-NS	71
FOB · FOB-NL	72
FXN · FXN-M	75

G	
GAB-N	64
H	
HCB	171
HCD	173
HCN	172
HDB-N · HDB-XL	170
HFA	174
HFC	175
HFD	176-177
HFM	180

I	
IFB	114
IFC NEW	112-113
IFS	92
IOC	110-111
IXC NEW	116-117

K	
KBP-V20 IP65 NEW	152
KDP	139
KFA-V20	154
KFB-TM	148
KFC-TM	150
KFD-V20	153
KFD-V40 IP68	155
KFN-TM	149
KFP-V20 IP65	152
KFP-V20 IP67	153
KFP-V30	154
KFP-V40 IP68 NEW	155
KFS-TM	148
KFU-V20	154
KFU-V30	154
KGP	136
KIB-TM	149
KIP-V20M	153
KXC-TM	151
KXP-V20 IP65	152

N	
NFB	133
NIB	132
NFN	134

P	
PBJ/PBS	39
PCB	28-29
PCD	26
PCJ NEW	32-33
PEJ/PES	40
PFB	27
PLJ/PLS	35
PNS/PNJ	36
PWS	38

R	
RIB	80
RPB	81

S	
SCD-4.0	186
SCD-4.0-PRO	187
SFB/SFB-H	78
SFE	118
SXC NEW	76-77

T	
TGC	11
TGD	11

U	
UFA	127
UFB	130
UFN	131
UIB	128
UID	129

V	
VHB	165

W	
WTB	73

Y	
YKV	137
YRO-01/-02/-03	192



KERN Quick-Finder

How I quickly find the product I am looking for

Ahead of each product group allows you to base the search for a certain target group on weighing data you need such as readout, weighing capacity and main features for each model.

And it's as simple as that – find the product you want in 2 steps:

1. Go to the product group index on page 3
2. Pick the appropriate product group and find the product you want using the Quick-Finder



The price leader in analytical balances, with internal or external adjustment

Features

- KERN ADJ: Automatic internal adjustment in the case of a change in temperature $\geq 2\text{ }^{\circ}\text{C}$ or timecontrolled every 3 h, guarantees high degree of accuracy and makes the balance independent of its location of use
- KERN ADB: Adjusting program CAL for quick setting of the balance accuracy using an external test weight at an additional price, see *Test Weights*
- Glass draught shield with 3 sliding doors for easy access to the items being weighed. standard
- **1** ADB 600-C3/ADJ 600-C3: Compact, space-saving carat balances with a readout of 0.001 ct and a weighing range of 600 ct. The high level of accuracy saves hard cash wherever you are weighing valuable precious stones
- Level indicator and levelling feet for precise levelling of the scale, fitted as standard, to give the most accurate weighing result
- Simple and convenient 6-key operation


Technical data

- Large backlit LCD display, Digit height 16 mm
- Dimensions weighing surface, stainless steel, $\varnothing 90\text{ mm}$
- Weighing space W×D×H
KERN ADB-C/ADJ-C: 170×160×110 mm
KERN ADB/ADJ: 170×160×205 mm
- Permissible ambient temperature 10 °C/30 °C

Accessories

- **2** Set for density determination of liquids and solids with density $\leq/\geq 1$, the density is indicated directly on the display, KERN YDB-03, **€ 380,-**
- **3** Ioniser to neutralise electrostatic charge, KERN YBI-01A, **€ 1000,-**
- **4** Gemstones plate, aluminium with practical spout, W×D×H 130×80×30 mm, KERN AEJ-A05, **€ 43,-**
- **5** Weighing table to absorb vibrations and oscillations, which would otherwise distort the weighing result, KERN YPS-03, **€ 1030,-**
- Minimum weight of sample, smallest weight to be weighed, depending on the required process accuracy, only in combination with a DAkkS calibration certificate, KERN 969-103, **€ 10,-**
- Equipment qualification: compliant qualification concept which includes the following validation services, Installation Qualification (IQ), Operating Qualification (OQ)
- Further details, plenty of further accessories and suitable printers see *Accessories*


STANDARD

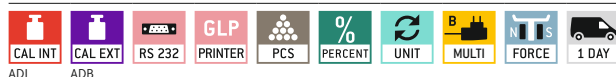

OPTION





Model	Weighing capacity [Max]	Readability [d]	Reproducibility	Linearity	Overall dimensions W×D×H mm	Net weight kg	Price excl. of VAT ex works €	Options DAKKS Calibr. Certificate
KERN	g	mg	mg	mg				DAKKS KERN €
ADB 100-4	120	0,1	0,2	± 0,4	230×310×330	4,4	770,-	963-101 191,-
ADB 200-4	220	0,1	0,2	± 0,4	230×310×330	4,4	850,-	963-101 191,-
ADB 600-C3	600 ct	0,001 ct	0,002 ct	± 0,004 ct	230×310×210	3,8	760,-	963-101 191,-
ADJ 100-4	120	0,1	0,2	± 0,4	230×310×330	5,0	910,-	963-101 191,-
ADJ 200-4	220	0,1	0,2	± 0,4	230×310×330	5,0	990,-	963-101 191,-
ADJ 600-C3	600 ct	0,001 ct	0,002 ct	± 0,004 ct	230×310×210	4,6	900,-	963-101 191,-

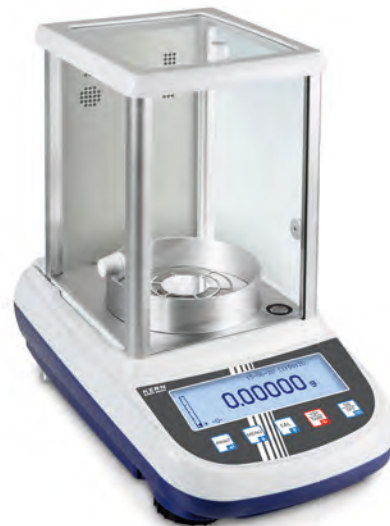

KERN ALJ 200-5DA with optional ioniser **I**, see *accessories*. High-precision semi-micro analytical balance. Thanks to its high level of precision, it is ideal for calibrating pipettes
 Note: To prevent evaporation we recommend economical capillary tubes (see standard 8655)

Analytical balances with a large weighing range, graphics display and user-friendly recipe weighing function – also as single-range semi-micro balance with unbelievably high resolution

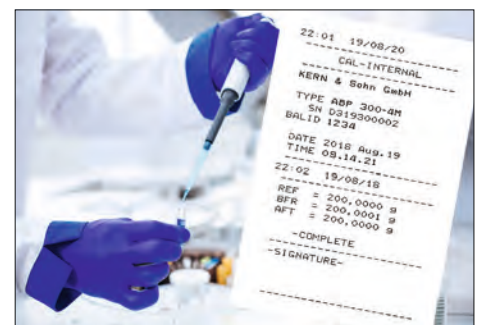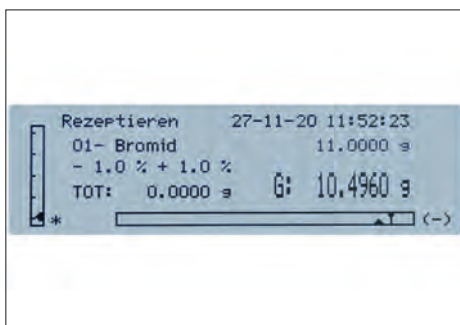

Convenient recipe-weighing: with the recipe database, up to 99 recipes can be stored, each with up to 20 recipe ingredients with name and target value

Clear printout with date and time. In addition, the components of the mixture are numbered automatically and printed out with the name and weight

GLP/ISO record keeping: professional, detailed GLP Protocol, so that the balance is completely compliant with the relevant standard requirements in accordance with ISO, GLP and GMP



Features

- ALJ 210-5A: Semi-micro model with just one weighing range with unbelievably high resolution, ideal where heavy items need to be weighed with the most accurate readout across the entire weighing range. **1** Particularly advantageous: the ioniser KERN ALJ-A03 for neutralising electrostatic charge is already fitted as standard
- Rapid and efficient operation thanks to the graphics display. Simple, clear user interface on the display in the following languages: DE, EN, FR, IT, ES, PT
- KERN ALJ: Automatic internal adjustment, guarantees high degree of accuracy and makes the balance independent of its location of use
- KERN ALS: Adjusting program CAL for quick setting of the balance accuracy using an external test weight at an additional price, see *Test Weights*
- Short stabilisation time: steady weight values within approx. 4 s (Models with [d] = 0,1 mg) 10 s; 6 s (Models with [d] = 0,01 | 0,1 mg) under laboratory conditions
- Weighing with tolerance range (checkweighing): Input of an upper/lower limit value. A visual signal assists with portioning, dispensing or grading
- Dosage aid: High stability mode and other filter settings can be selected

- Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display
- Ergonomically optimised keypad for left and righthanded users
- Large glass draught shield with 3 sliding doors for easy access to the items being weighed. Weighing space W×D×H 160×170×225 mm
- Compact size, practical for small spaces
- Protective working cover included with delivery

Technical data

- Backlit LCD graphic display, digit height 15 mm
- Dimensions weighing surface, stainless steel, ø 80 mm
- Overall dimensions W×D×H 210×340×330 mm
- Net weight approx. 7 kg
- Permissible ambient temperature 5 °C/35 °C

Accessories

- Protective working cover, scope of delivery 5 items, KERN ALJ-A01S05, **€ 44,-**
- Protective dust cover, KERN ABS-A08, **€ 32,-**
- 1** Evaporation trap, minimises faults through evaporation when using pipettes for small volumes of 10 µl to 10 ml, KERN ALJ-A02, **€ 720,-**
- 2** Draught shield rear panel with integrated ioniser to neutralise electrostatic charge. Particularly convenient handling as you no longer need a separate device. Simply enable the ionizer fan at the push of a button.. Is fitted in place of the existing glass rear panel of the draught shield. Please order at the time you order your balance, the scope of delivery is the rear panel, ionizer, Universal plug-in power supply. Factory Option, KERN ALJ-A03, **€ 620,-**
- 3** Set for density determination of liquids and solids with density ≤/≥ 1, the density is indicated directly on the display, KERN YDB-03, **€ 380,-**
- Weighing table to absorb vibrations and oscillations, which would otherwise distort the weighing result, KERN YPS-03, **€ 1030,-**
- Equipment qualification: compliant qualification concept which includes the following validation services, Installation Qualification (IQ), Operating Qualification (OQ)
- Further details, plenty of further accessories and suitable printers see *Accessories*

STANDARD

CAL INT	CAL EXT	RS 232	GLP	INTERN	PCS	RECIPE	PERCENT	UNIT	TOL	MOVE	MULTI	FORCE	1 DAY
---------	---------	--------	-----	--------	-----	--------	---------	------	-----	------	-------	-------	-------

ALJ ALS

OPTION FACTORY

DAkKS	M
+3 DAYS	ALJ-M

Model	Weighing capacity [Max] g	Readability [d] mg	Verification value [e] mg	Minimal load [Min] mg	Reproducibility mg	Linearity mg	Price excl. of VAT ex works €	Verification M T KERN €	Options DAkKS Calibr. Certificate DAkKS KERN €
ALS 160-4A	160	0,1	-	-	0,1	± 0,3	1330,-	-	963-101 191,-
ALS 250-4A	250	0,1	-	-	0,1	± 0,3	1380,-	-	963-101 191,-
ALJ 210-5A	210	0,01	-	-	0,05	± 0,1	3300,-	-	963-101 191,-
ALJ 160-4A	160	0,1	-	-	0,1	± 0,3	1420,-	-	963-101 191,-
ALJ 250-4A	250	0,1	-	-	0,1	± 0,3	1450,-	-	963-101 191,-
ALJ 310-4A	310	0,1	-	-	0,1	± 0,3	1500,-	-	963-101 191,-
ALJ 500-4A	510	0,1	-	-	0,2	± 0,4	2250,-	-	963-101 191,-

Multi-range balance, with increasing load it switches automatically to the next largest weighing range [Max] and readout [d] and when the load is fully removed, the balance switches back to the lower range


ALJ 200-5DA	82 220	0,01 0,1	-	-	0,04 0,1	± 0,1 0,2	2250,-	-	963-101 191,-
--------------------	----------	------------	---	---	------------	-------------	---------------	---	---------------

Note: For devices that require verification (conformity assessment according to NAWI 2014/31/EU), please include the verification when placing your order.

The initial verification is not possible after delivery. Please inform the full address of the location of use for the initial verification.

ALJ 160-4AM	160	0,1	1	10	0,1	± 0,3	1460,-	965-201 150,-	963-101 191,-
ALJ 250-4AM	250	0,1	1	10	0,1	± 0,3	1490,-	965-201 150,-	963-101 191,-

1 ONLY WHILE STOCKS LAST

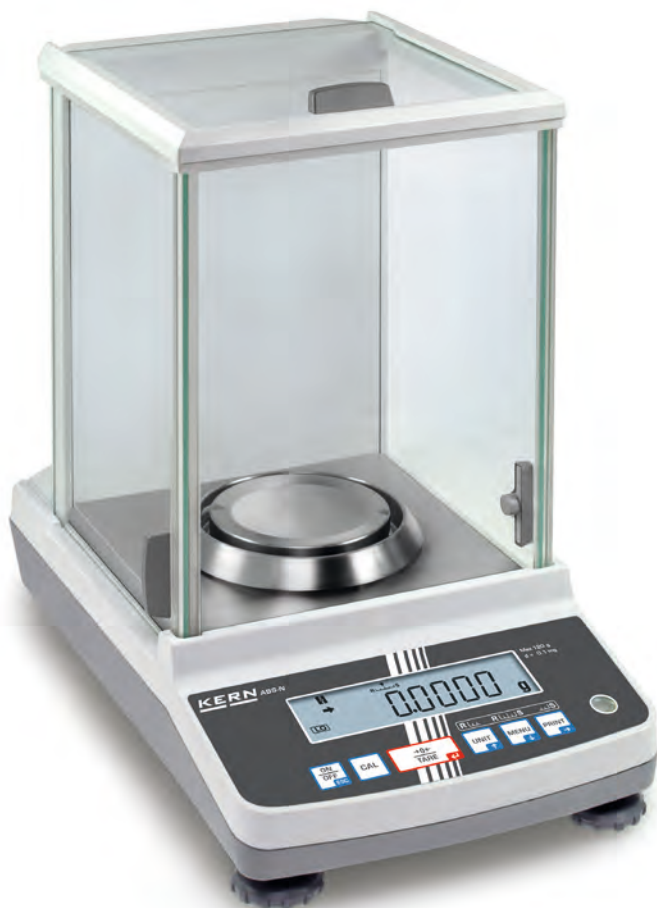

USB data interfaces and RS-232 for transferring weighing data to the PC, tablet, printer or USB

The bestseller in analytical balances, with high-quality single-cell weighing system, verification optional

Features

- KERN ACJ: Automatic internal adjustment in the case of a change in temperature \geq or timecontrolled every , guarantees high degree of accuracy and makes the balance independent of its location of use
- KERN ACS: Adjusting program CAL for quick setting of the balance accuracy using an external test weight at an additional price, see *Test Weights*
- Dosage aid: High stability mode and other filter settings can be selected
- Simple recipe weighing and documenting with a combined tare/print function. In addition, the ingredients for the recipe are numbered automatically and printed out with their corresponding number and nominal weight
- Automatic data output to the PC/printer each time the balance is steady
- Identification number: 4 digits, printed on calibration protocol freely programmable
- Protective working cover included with delivery



Technical data

- Large LCD display, digit height 14 mm
- Dimensions weighing surface, stainless steel, \varnothing 91 mm
- Weighing space W×D×H 174×162×227 mm
- Overall dimensions W×D×H 213×333×338 mm
- Net weight approx. 7 kg
- Permissible ambient temperature 10 °C/30 °C

Accessories

- Protective working cover, scope of delivery 5 items, KERN ACS-A02S05, € 44,-
- **1** Set for density determination of liquids and solids with density \leq/\geq 1, the density is indicated directly on the display, KERN YDB-03, € 380,-
- **2** Ioniser to neutralise electrostatic charge, KERN YBI-01A, € 1000,-
- **3** Weighing table to absorb vibrations and oscillations, which would otherwise distort the weighing result, KERN YPS-03, € 1030,-
- Minimum weight of sample, smallest weight to be weighed, depending on the required process accuracy, only in combination with a DAkkS calibration certificate, KERN 969-103, € 10,-
- Equipment qualification: compliant qualification concept which includes the following validation services, Installation Qualification (IQ), Operating Qualification (OQ)
- Further details, plenty of further accessories and suitable printers see *Accessories*



- 4** Single-cell advanced technology:
- Fully automatic manufactured weighing cell from one piece of material
 - Stable temperature behaviour
 - Short stabilisation time: steady weight values within approx. 3 s under laboratory conditions
 - Shock proof construction
 - High corner load performance

STANDARD

ACJ	ACS												

OPTION

+3 DAYS	+3 DAYS

FACTORY

ACJ

Model	Weighing capacity [Max]	Readability [d]	Verification value [e]	Minimal load [Min]	Reproducibility	Linearity	Price excl. of VAT ex works €	Options			
								Verification	€	DAkkS	Certificate
KERN	g	mg	mg	mg	mg	mg	€	MET	€	DAkkS	€
ACS 80-4	82	0,1	-	-	0,2	± 0,3	1230,-	-	-	963-101	191,-
ACS 100-4	120	0,1	-	-	0,2	± 0,3	1330,-	-	-	963-101	191,-
ACS 200-4	220	0,1	-	-	0,2	± 0,3	1430,-	-	-	963-101	191,-
ACS 300-4	320	0,1	-	-	0,2	± 0,3	1630,-	-	-	963-101	191,-

Note: For devices that require verification (conformity assessment according to NAWI 2014/31/EU), please include the verification when placing your order. The initial verification is not possible after delivery. Please inform the full address of the location of use for the initial verification.

Model	Weighing capacity	Readability	Verification	Minimal load	Reproducibility	Linearity	Price	Verification	€	DAkkS	Certificate
ACJ 80-4M	82	0,1	1	10	0,2	± 0,3	1370,-	965-201	150,-	963-101	191,-
ACJ 100-4M	120	0,1	1	10	0,2	± 0,3	1430,-	965-201	150,-	963-101	191,-
ACJ 200-4M	220	0,1	1	10	0,2	± 0,3	1530,-	965-201	150,-	963-101	191,-
ACJ 300-4M	320	0,1	1	10	0,2	± 0,3	1730,-	965-201	150,-	963-101	191,-



4

The premium model with single-cell weighing system

Features

- Automatic internal adjustment in the case of a change in temperature $\geq 0,5\text{ }^{\circ}\text{C}$ or timecontrolled every 4 h, guarantees high degree of accuracy and makes the balance independent of its location of use
- Dosage aid: High stability mode and other filter settings can be selected
- Simple recipe weighing and documenting with a combined tare/print function. In addition, the ingredients for the recipe are numbered automatically and printed out with their corresponding number and nominal weight
- Identification number: 4 digits, printed on calibration protocol freely programmable
- Printout of a GLP-compliant calibration report conveniently at the touch of a button
- Automatic data output to the PC/printer each time the balance is steady

- Large glass draught shield with 3 sliding doors for easy access to the items being weighed.
- Protective working cover included with delivery

Technical data

- Large LCD display, digit height 14 mm
- Dimensions weighing surface \varnothing 80 mm, stainless steel
- Weighing space WxDxH 168x172x223 mm
- Overall dimensions WxDxH 217x356x338 mm
- Net weight approx. 7 kg
- Permissible ambient temperature $10\text{ }^{\circ}\text{C}/30\text{ }^{\circ}\text{C}$


Single-cell advanced technology:

- Fully automatic manufactured weighing cell from one piece of material
- Stable temperature behaviour
- Short stabilisation time: steady weight values within approx. 4 s (models with [d] = 0,1 mg), approx. 10 s (models with [d] = 0,01 mg) under laboratory conditions
- Shock proof construction
- High corner load performance

Accessories

- Protective working cover, scope of delivery 5 items, KERN ABT-A02S05, **€ 44,-**
- **1** Set for density determination of liquids and solids with density $\leq/\geq 1$, the density is indicated directly on the display, KERN YDB-03, **€ 380,-**
- **2** Ioniser to neutralise electrostatic charge, KERN YBI-01A, **€ 1000,-**
- **3** Weighing table to absorb vibrations and oscillations, which would otherwise distort the weighing result, KERN YPS-03, **€ 1030,-**
- Minimum weight of sample, smallest weight to be weighed, depending on the required process accuracy, only in combination with a DAkkS calibration certificate, KERN 969-103, **€ 10,-**
- Equipment qualification: compliant qualification concept which includes the following validation services, Installation Qualification (IQ), Operating Qualification (OQ)
- Further details, plenty of further accessories and suitable printers see *Accessories*

STANDARD

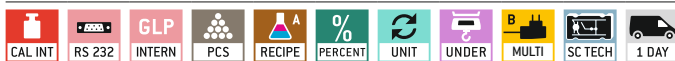

OPTION



FACTORY



Model	Weighing capacity [Max] g	Readability [d] mg	Verification value [e] mg	Minimal load [Min] mg	Reproducibility mg	Linearity mg	Price excl. of VAT ex works €	Options			
								Verification MID KERN	€	DAkkS Calibr. Certificate DAkkS KERN	€
ABT 100-5NM	101	0,01	1	1	0,05	$\pm 0,15$	2500,-	965-201	150,-	963-101	191,-
ABT 120-4NM	120	0,1	1	10	0,1	$\pm 0,2$	1900,-	965-201	150,-	963-101	191,-
ABT 220-4NM	220	0,1	1	10	0,1	$\pm 0,2$	2300,-	965-201	150,-	963-101	191,-
ABT 320-4NM	320	0,1	1	10	0,1	$\pm 0,3$	2600,-	965-201	150,-	963-101	191,-
Multi-range balance, with increasing load it switches automatically to the next largest weighing range [Max] and readout [d] and when the load is fully removed, the balance switches back to the lower range											
ABT 120-5DNM	42 120	0,01 0,1	1 1	1	0,02 0,1	$\pm 0,05 0,2$	2550,-	965-201	150,-	963-101	191,-
ABT 220-5DNM	82 220	0,01 0,1	1 1	1	0,05 0,1	$\pm 0,1 0,2$	2650,-	965-201	150,-	963-101	191,-

Note: For devices that require verification (conformity assessment according to NAWI 2014/31/EU), please include the verification when placing your order. The initial verification is not possible after delivery. Please inform the full address of the location of use for the initial verification.



2 KERN ABP 100-5DM
with optional ioniser

Premium analytical balance with the latest Single-Cell Generation for extremely rapid, stable weighing results - now also as a version with automatic sliding doors



1 The internal draught shield minimises the effect of currents of air in the weighing chamber and therefore significantly improves the stabilisation time and repeatability

2 Draught shield rear panel with integrated ioniser, which can be fitted in place of the existing glass rear panel of the draught shield. Extremely fast ionization process, thanks to the latest generation of KERN ionization technology to neutralise electrostatic charge

GLP/ISO record keeping: professional, detailed GLP Protocol, so that the balance is completely compliant with the relevant standard requirements in accordance with ISO, GLP and GMP

Features

- This new generation of analytical balances combines the highest level of precision with large weighing ranges. Thanks to the new Single-Cell Generation, the weighing result is displayed in a fraction of the time with comparable models. Together with the intuitively structured menu, this means that you can work efficiently and rapidly
- Large glass draught shield with 3 sliding doors for easy access to the items being weighed
- KERN ABP-A: with 3 automatic sliding doors, which can be opened and closed using sensors, and integrated ioniser as standard. Thanks to the Memory Function, the balance stores how far the sliding doors have been opened. This prevents contamination and accelerates processes. Touching the sliding doors “gently” activates the push function and the doors open and shut automatically. The adjustable internal draught shield guarantees maximum stability of the weight readings
- Navigation pad for super quick navigating through the menus
- Automatic internal adjustment in the case of a change in temperature $\geq 1\text{ }^{\circ}\text{C}$ or timecontrolled every 4 h, guarantees high degree of accuracy and makes the balance independent of its location of use
- The minimum weight of sample can be manually stored in the device or automatically calculated. For weighings below this value, the balance issues a warning message
- Dosage aid: High stability mode and other filter settings can be selected
- Simple recipe weighing and documenting with a combined tare/print function. In addition, the ingredients for the recipe are numbered automatically and printed out with their corresponding number and nominal weight
- Individual user settings can be stored for up to 10 users: user name/number (can be printed out or added to the record for each process),

- password, menu language, user profiles, accessing user settings via barcode, additional guest mode for users who are not logged in, authorizations, eg. B. balance adjustment, changing settings or conditioning or modification of a recipe only by the authorized person & performing the formulation by the user
- RS 232 data interfaces and USB (device) for transferring weighing data and USB (host) to connect a USB keyboard for easy capture of item numbers, recipes, for easier navigation within the menu etc.
- U.S. FDA 21 Part 11: assists you in data integrity in accordance with U.S. FDA 21 Part 11 (for example weighing result, sample ID, user name, scales ID, ...)
- Menu languages DE, EN
- Automatic data output to the PC/printer each time the balance is steady for models with [d] = 0,0001 g
- Multi-function weighing plate included with delivery, minimises the effect of currents of air in the weighing chamber and therefore significantly improves the stabilisation time and repeatability. In addition samples, sample paper, PCR containers, micro centrifuge tubes and many other items which protrude can be easily fixed in place and weighed easily
- ABP 200-5M: Erlenmeyer flask holder included with delivery
- Protective working cover included with delivery

Technical data

- Luminescent OLED display, digit height 12 mm bright with high contrast, for easy reading of the weight, even in poor lighting conditions
- Dimensions weighing surface \varnothing 91 mm, stainless steel
- Weighing space WxDxH 166x156x220 mm
- Overall dimensions WxDxH 220x370x350 mm
- Net weight approx. 8 kg
- Permissible ambient temperature 10 °C/30 °C


Accessories

- Protective working cover, scope of delivery 5 items, KERN YBA-A06S05, € 44,-
- Set for density determination of liquids and solids with density $\leq/\geq 1$, the density is indicated directly on the display, KERN YDB-03, € 380,-
- **1** Internal draught shield made of glass, KERN ABP-A02, € 165,-
- **2** Draught shield rear panel with integrated ionizer to neutralise electrostatic charge. Is fitted in place of the existing glass rear panel of the draught shield. Suitable for all models in the range, please order at the time you order your balance, the scope of delivery is the rear panel, ionizer, Universal plug-in power supply. Factory Option, KERN ABP-A01, € 1990,-
- USB barcode scanner, hand-held version, dimensions WxDxH 152x84x63 mm, KERN PET-A09, € 300,-
- Minimum weight of sample, smallest weight to be weighed, depending on the required process accuracy, only in combination with a DAkKS calibration certificate, KERN 969-103, € 10,-
- Equipment qualification: compliant qualification concept which includes the following validation services, Installation Qualification (IQ), Operating Qualification (OQ)
- Further details, plenty of further accessories and suitable printers see *Accessories*

Single-cell advanced technology:

- Fully automatic manufactured weighing cell from one piece of material
- Stable temperature behaviour
- Short stabilisation time: steady weight values within approx. 2 s (models with [d] = 0,1 mg), approx. 8 s (models with [d] = 0,01 mg) under laboratory conditions
- Shock proof construction
- High corner load performance

STANDARD

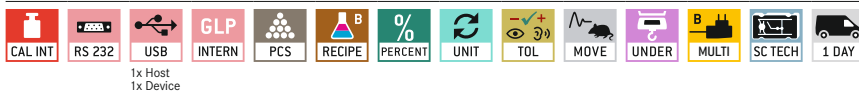

OPTION



FACTORY



Model	Weighing capacity [Max] g	Readability [d] mg	Verification value [e] mg	Minimal load [Min] mg	Reproducibility mg	Linearity mg	Price excl. of VAT ex works €	Options			
								Verification KERN	€	DAkKS Calibr. Certificate DAKKS KERN	€
ABP 100-5M	135	0,01	1	1	0,05	$\pm 0,1$	4200,-	965-201	150,-	963-101	191,-
ABP 200-5M	220	0,01	1	1	0,05	$\pm 0,1$	6000,-	965-201	150,-	963-101	191,-
ABP 100-4M	120	0,1	1	10	0,1	$\pm 0,2$	2200,-	965-201	150,-	963-101	191,-
ABP 200-4M	220	0,1	1	10	0,1	$\pm 0,2$	2700,-	965-201	150,-	963-101	191,-
ABP 300-4M	320	0,1	1	10	0,2	$\pm 0,3$	3000,-	965-201	150,-	963-101	191,-
Multi-division balance, with increasing or decreasing load, it switches automatically to the next largest or smallest weighing range [Max] and readout [d].											
ABP 100-5DM	52 120	0,01 0,1	1 1	1	0,02 0,1	$\pm 0,05$ 0,2	3200,-	965-201	150,-	963-101	191,-
ABP 200-5DM	102 220	0,01 0,1	1 1	1	0,05 0,1	$\pm 0,1$ 0,2	3400,-	965-201	150,-	963-101	191,-

Note: For devices that require verification (conformity assessment according to NAWI 2014/31/EU), please include the verification when placing your order. The initial verification is not possible after delivery. Please inform the full address of the location of use for the initial verification.



Automatic sliding doors, which can be opened and closed using sensors



Ioniser (included) can be activated at the touch of a button

4



Opening possible alternatively via button



KERN ABP-A with integrated ioniser and internal draft shield



Erlenmeyer piston holder, included as standard with ABP 200-5M and ABP 200-5AM




Internal draft shield - minimises the effect of air currents in the weighing chamber, only for 0.01 mg models of the series ABP-A



Multi-function weighing plate, Note: included with delivery for models with [d] = 0,01 mg

STANDARD

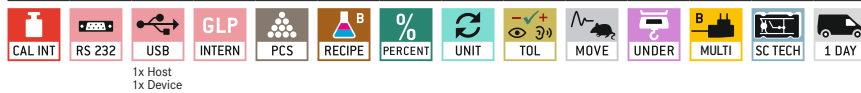

OPTION



FACTORY



Model	Weighing capacity [Max]	Readability	Verification value	Minimal load	Reproducibility	Linearity	Price excl. of VAT ex works €	Options			
								Verification	DAkKS Calibr.	Certificate	
KERN	g	[d] mg	[e] mg	[Min] mg	mg	mg	€	€	KERN	€	
ABP 100-5AM	135	0,01	1	1	0,05	± 0,1	9200,-	965-201	150,-	963-101	191,-
ABP 200-5AM	220	0,01	1	1	0,05	± 0,1	11200,-	965-201	150,-	963-101	191,-
ABP 200-4AM	220	0,1	1	10	0,1	± 0,2	5900,-	965-201	150,-	963-101	191,-
ABP 300-4AM	320	0,1	1	10	0,2	± 0,3	6400,-	965-201	150,-	963-101	191,-

Multi-division balance, with increasing or decreasing load, it switches automatically to the next largest or smallest weighing range [Max] and readout [d].

ABP 100-5DAM	52 120	0,01 0,1	1 1	1	0,02 0,1	± 0,05 0,2	6800,-	965-201	150,-	963-101	191,-
ABP 200-5DAM	102 220	0,01 0,1	1 1	1	0,05 0,1	± 0,1 0,2	9000,-	965-201	150,-	963-101	191,-

Note: For devices that require verification (conformity assessment according to NAWI 2014/31/EU), please include the verification when placing your order. The initial verification is not possible after delivery. Please inform the full address of the location of use for the initial verification.

KERN & SOHN GmbH

Balances, Test Weights, Microscopes,
DAkkS Calibration Laboratory
Ziegelei 1
72336 Balingen
Germany
Tel. +49 7433 9933-0
info@kern-sohn.com
www.kern-sohn.com

**The oldest Precision Balance
Factory in Germany**



**Discover the multifaceted World of Balances, Microscopes and Measuring Technology
from KERN online: www.kern-sohn.com**



- Full KERN & SAUTER Product Range
- Convenient 24/7 Ordering
- Selection of more than 5,000 Items across Weighing and Measuring Technology, Optical Instruments as well as Accessories and Services
- Extensive Information and useful Download Options
- Technical Product Data Sheets
- Operating Instructions
- Descriptive Image and Video Material
- Useful KERN Services
- Technical Glossary
- KERN Dealer Portal
- Practical Filter and Search Functions



Follow us on our social media channels too



Printed in Germany by KERN & SOHN GmbH
z-cb-en-kr-20241

